

## DIRECTIONS FOR USE

### EDGEFILE®X3 Heat Treated Fire-Wire™

#### Intended Purpose

Endodontic files and reamers are single use surgical instruments used for performing root canal treatment to mechanically shape and prepare the root canals during endodontic therapy or to remove the root canal obturating material when performing retreatment. The device is intended to be used sterile and single use only.

#### COMPOSITION

The instrument is made of a nickel-titanium blade, handle, the stop, and the color-coded band.

#### Contraindications

- Mechanically driven endodontic instruments should not be used in cases with very severe and sudden curvatures.
- This product contains nickel and should not be used for individuals with known allergic sensitivity to this metal.

#### Warnings

- A rubber dam system should be used.
- Endodontic files are sharp, and caution should be used if touching the blade directly.
- The endodontic files must be sterilized before patient use.
- Endodontic files are for single use only, in order to avoid file separation.
- This device is intended for single patient use only.
- Used files should be disposed of in a Biohazard Sharps container in accordance with local regulations.
- After use, this product may be a potential biohazard. Handle and dispose of in accordance with accepted medical practice and applicable local, state, and federal laws and regulations.
- Do not use if package is damaged.
- The product and the packaging have not been designed or tested for reuse. The ability to effectively clean and re-sterilize this single use device and subsequent reuse may adversely affect the clinical performance, safety and/or sterility of the device.
- If the packaging is damaged, please dispose of the damaged product and utilize an undamaged product instead as the former may be contaminated.

#### Precautions for Use

As with all products, use carefully until you become proficient with use. Always determine working length using radiographs and/or apex locator to properly use endodontic files. Important points to remember:

1. Use only in an electric motor and hand piece designed for endodontic (rotary/reciprocating) files.
2. Straight-line access is imperative for proper file use and endodontic treatment.
3. Do not force the files down canals, use minimal apical pressure.
4. Clean the flutes frequently and at least after removing the files from the canal.
5. Irrigate and lubricate the canal frequently throughout the procedure.

6. Take each file to length only one time and for no more than one second.
7. In apical areas and curved canals exercise caution.
8. Endodontic files are single patient use devices.
9. Once file is used do not reuse. If file is reused and used on a different patient infection can be introduced. Performance of the file can also be reduced.
10. When instrumenting the canal, do not over enlarge the coronal portion of the canal.
11. Too large a file taken to length increases the risk of canal transportation and file separation.
12. Do not exceed the handpiece recommended maximum torque or speed. Exceeding settings may cause the device to fail.
13. Endodontic files undergo our proprietary Annealed Heat Treatment (AHT) forming our branded Fire-Wire™ NiTi which increases cyclic fatigue resistance and torque strength. With this proprietary processing, the files may be slightly curved. This is not a manufacturing defect. While the file can be easily straightened with your fingers, it is not necessary as once they are inside the canal, endodontic files will follow and conform to the natural canal anatomy and curvatures.
14. Do not use after the expiration date on the label.

#### Adverse Reactions

- Device fracture/breakage
- Infection – Do not use if package is damaged or open, due to risk of infection occurring.
- Complications usually associated with endodontic procedures including:
  - Pain
  - Instrument fracture/breakage
  - Soft tissue damage/bleeding

#### Safe Unwinding

As a safety feature the files are designed to unwind. They may be used until the files unwind backwards.

## INSTRUCTIONS FOR USE

#### Sterilisation

No sterilisation steps are needed for product provided sterile.

#### Straight-Line Access and Glide Path Formation

- Prepare straight-line access to all canal orifice.
- With lubrication in the canal form a glide path with a size #10 and #15 hand files or mechanical glide path files 2/3 down the length of the canal.

#### Beginning Canal Shaping and Cleaning: N1 and N2 Files

- With lubricant in the canal and with light apical pressure take the N1 (17/06) into the canal and follow the glide path using an in-and-out motion while laterally brushing the dentin on the outstroke to enhance the straight-line access of the canal.
- Continue shaping with the N1 until resistance is met or 2/3 down the canal is reached.
- Then use the N2 (17/04), in the same way until resistance is met or 2/3 down the canal is reached.

DFU-SEFX3-EU Rev A

- Switch between the N1 and N2 following the glide path using the same in-and-outas described for both files until 2/3 down the canal is reached.
- Now that the coronal 2/3 of the canal is shaped, form a glide path with the size #10 and #15 hand files or mechanical glide path files into the apical 1/3.
- Establish working length with radiographs and/or an apex locator. Then confirm patency by taking the #10 hand file 1 mm past the working length.
- Then, using the same motion as before, switch between the N1 (17/06) and N2 (17/04) until N2 reaches the working length.
- If a larger coronal shape is desired, use the NX (25/12) at any time after the coronal 2/3 is shaped. Completing Canal Shaping and Cleaning: C1, C2, C3, C4 Files
- With lubricant in the canal and with light apical pressure complete canal shaping and cleaning by taking the C1 (20/06) down the canal until the working length is reached.
- Apically gauge the foramen at the working length with a #20 hand file. If the #20 hand file is snug at the working length, the canal is shaped and ready to obturate.
- If the #20 hand file is loose, take the C2 (25/06) to the working length, then gauge with a #25 hand file. When necessary, the C3 (30/06) or C4 (40/06) may need to be used.

### Motor Settings

The **EDGEFILE®X3** can be used in a clockwise reciprocating motor but not in the WaveOne reciprocating motor, using the WaveOne setting, which moves in the counterclockwise direction. The **EDGEFILE®X1** is designed specifically for use in only the WaveOne reciprocating motor and setting.

Use the same hand piece with the same speed and torque settings you are currently using with your ProTaper® or ProTaper Next® rotary system. Or if you wish, you can use for all **EDGEFILE®X3** rotary files the following speed and torque settings for all files.

- Speed: 300-500 rpm
- Torque: 2.94 N-cm/300 g-cm

### Disinfecting

- After each canal is fully shaped, rinse the canals for 1 minute with 17% Liquid EDTA to remove the canal Smear Layer.
- Rinse the canals for 5 minutes with 5% NaOCl to remove debris and bacteria.
- Rinse the canals for 1 minute with 17% Liquid EDTA to rinse out the 5% NaOCl.
- Rinse the canals for 5 minutes with 2% chlorohexidine or EDTA to kill bacteria.

### Obturation of Canal Systems

- When using thermal carrier system use size verifiers to determine the proper sized carrier.
- When using a master gutta percha cone that matches the largest file taken to length, remember sometimes you may need to drop down in cone tip size if the corresponding gutta percha to your final rotary file does not go to length.





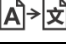


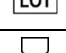











### Storage Conditions

- Store at room temperature of 10°C~37.8°C, away from any sunlight. Recommended File Disposal Place used files in Biohazard Sharps container.

### Reporting of Incidents to manufacturer and competent authorities

- In case any patient/user faces a serious incident, the entirety of the incident will be reported to the following:
  - The manufacturer of the device: US ENDODONTICS
  - The competent authority of the country where the user/patient resides

### Symbol Table

Symbol	Meaning (Standard, if Applicable)
	Manufacturer: Indicates the medical device manufacturer (ISO 15223-1)
	Authorized Representative: Indicates the AR in the EU
	Importer: Indicates the entity importing the medical device into the locale (ISO 15223-1)
	Conformité Européenne. EU mandatory conformity marking.
	Translation: Indicates that the original information has been translated and replaced (ISO 15223-1)
	Medical Device: Indicates the item is a medical device (ISO 15223-1)
	Catalogue number: indicates the Medical Device SKU (ISO 15223-1)
	Batch Code: Manufacturer's batch code so batch or lot can be identified (ISO 15223-1)
	Indicates the date after which the medical device should not be used (ISO 15223-1)
	Do not reuse: Indicates a medical device that is intended for one single use only (ISO 15223-1)
	Do not resterilise: Indicates medical devices that is not to be resterilised (ISO 15223-1)
	Medical device sterilized using irradiation and packaged with a single outer sterile barrier system (ISO 15223-1)
	Medical device should not be used if package is damaged and consult instructions for use (ISO 15223-1)
	Consult IFU: consult the Instructions For Use, and eIFU website listed (ISO 15223-1)
	Caution is necessary when operating device. Align cautions (ISO 15223-1)
	Temperature limit to which the medical device can be safely exposed (ISO 15223-1)
	Unique Device Identifier: Indicates a carrier that contains UDI information (ISO 15223-1)
	Date and Country of Manufacture: To identify the country of manufacture of products next to date of manufacture (ISO 15223-1)
	Prescription Use Only: Caution: Federal law restricts this device to sale by or on the order of a dentist (21CFR 801.109)

DFU-SEFX3-EU Rev A